



How to Use a Clave Multidose Vial Adapter

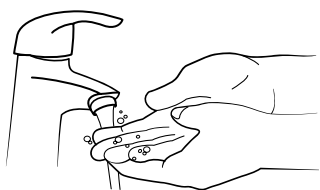


1.

Wash your hands well with soap and water, then dry them completely. Put on gloves if necessary.

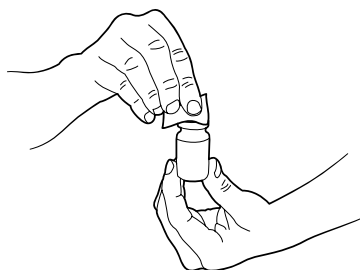


2.

Remove the plastic cap from the medication vial and use an alcohol wipe to clean the rubber stopper.

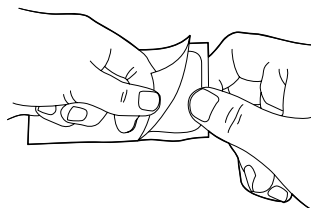


3.

Remove the multidose vial adapter from the package by holding the large rubber end. Remove the plastic tubing from the spike end and discard. Do not touch the now exposed spike end of the vial adapter.

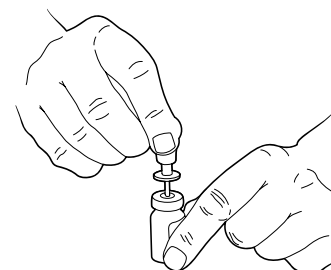

4.


Using the large rubber end to handle the adapter, connect the adapter to the vial by pushing the spike end straight down onto the top of the vial and in the center of the rubber stopper. Do not insert the spike at an angle. The spike will pierce through the rubber stopper. While pushing down on the vial adapter, turn $\frac{1}{2}$ rotation to create a tight seal on the medication vial's rubber stopper.

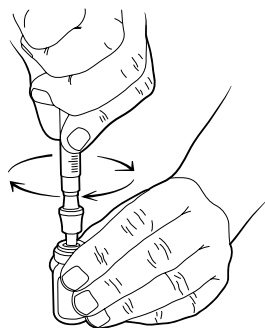

5.


Using a syringe without the needle attached, twist the syringe *clockwise* onto the large rubber end of the vial adapter to attach the syringe.

DO NOT use a syringe with a needle attached.

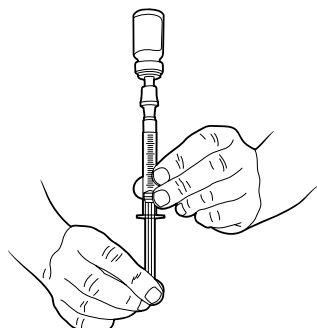


6.

Turn the vial, adapter, and syringe upside down to withdraw the prescribed amount for one injection.

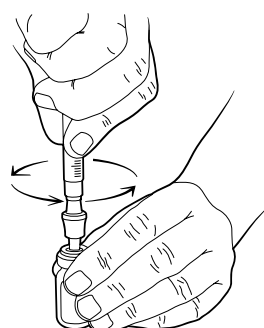

7.

Remove the syringe using a *counterclockwise* twisting motion. Confirm you have the right medication volume in the syringe. If not, simply reapply the syringe to the vial adapter and re-pull the medication to the right volume.

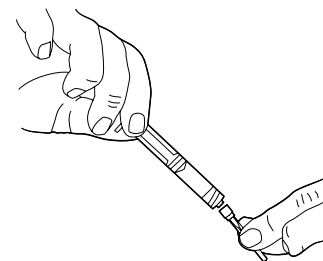

8.

Attach a capped needle onto the syringe using a clockwise twisting motion. Follow the injection technique described by your healthcare provider. Store the vial with the vial adapter attached.

DO NOT remove the vial adapter or reuse on another vial.